

STANDARD OPERATING PROCEDURES
SOP.UCF.036 - BD FACSYMPHONY A5 SE
STARTUP CLEANING PROCEDURE


Are you starting up?

(First user of the day on Weekends, Holidays and outside weekly working hours)

Turn on the Computer, the Cytometer and the FFSS. Login on Windows and turn ON all the lasers on the BD Coherent Connection 4.

Login on BD FACSDiva.

Follow the **startup cleaning procedure**:

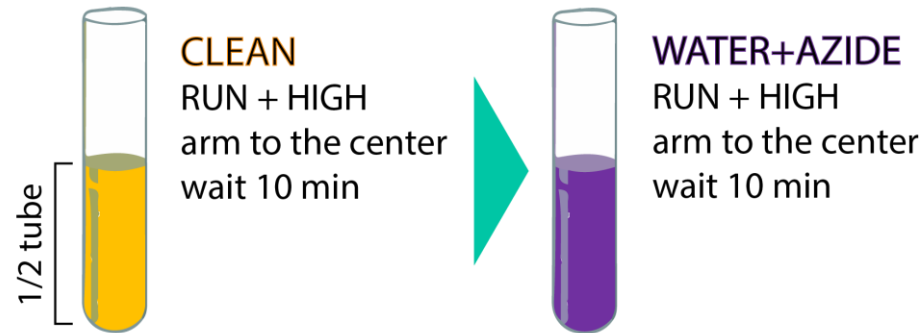

Perform the Prime Check with CS&T beads.
Run the Quality Control with CS&T beads.

You are ready to run your samples!

Thank you!